

AnimalTA user manual

v 2.3



Index

I.	Installation.....	4
II.	Generalities and project management	4
1.	Language	4
2.	Information panel	4
3.	Create a new project.....	4
4.	Open a project.....	4
5.	Save a project.....	4
III.	Gestion of videos.....	4
1.	Add videos.....	4
2.	Import videos from another project	5
3.	Suppress video	5
4.	Duplicate video.....	5
5.	Concatenate videos.....	5
IV.	Main project table	5
1.	Video identity	6
2.	Frame rate	6
3.	Cropping	7
4.	Stabilization	8
5.	Background.....	8
6.	Definition of the arenas	9
7.	Definition of the scale	10
V.	Video reader.....	11
VI.	Tracking preparation	12
	Original video	12
	Gray image	12
	Background.....	12
	Threshold.....	13
	Erosion/Dilation	13
	Area selection.....	14
	Distance threshold	14
	Number of individuals per arena	15
VII.	Tracking process.....	16
VIII.	Visualization and correction of the trackings.....	16

1.	Visualization	16
2.	Set the identities of the targets	17
3.	Corrections	17
	Find missing values.....	17
	Punctual error	17
	Interpolation	17
	Redo part of the tracking	18
	Invert coordinates	18
	Delete existing values.....	18
	Validate corrections	18
IX.	Prepare analyses	19
1.	Smoothing filter.....	19
2.	Perspective correction	19
3.	Analyses of individual trajectories	20
4.	Spatial repartition	21
	Elements manipulation	21
	Types of elements	21
	Expend elements to other arenas	22
5.	Interindividual distance.....	23
6.	Exploration	23
	Individual centered	23
	Meshes	23
X.	Run analyses.....	24
	Results per individual	24
	Results of interindividual distances.....	26
	Spatial.....	26
	Detailed data	27
XI.	Export a video.....	27

I. Installation

AnimalTA is a program developed in Python for usage under the Windows operating system. The AnimalTA installer can be downloaded from: <http://vchiara.eu/index.php/animalta> or <https://github.com/VioletteChiara/AnimalTA/releases>.

Download the file AnimalTA.installer.exe, launch it and follow the installer instruction.

II. Generalities and project management

1. Language

A drop menu is available on the top right corner of the AnimalTA window to change the displayed language. Note that changing the language will automatically close AnimalTA. Just reopen it to have the new set language.

2. Information panel

An information panel is always visible in the top right corner of the application. This panel aims to provide the user with information about the different components of the window simply by putting the mouse over them.

3. Create a new project

The notion of project is at the center of AnimalTA functioning. A project may be seen as a collection of videos, each associated with various information. By pressing the “New project” button, you will be able to choose a location on your computer where you want the files linked with your project to be located. Two things will be created in the indicated location: a *.ata* file which is your project’s file and a folder with the name of your project that will later contain all the outputs from AnimalTA.

4. Open a project

If you already have created a project before, you can reopen it using the “Open project button”. Then just select the *.ata* project file you want.

5. Save a project

An opened project can be saved using the “Save” Button. The last saved version of the project will be replaced by the current one.

III. Gestion of videos

1. Add videos

Use the “Add video(s)” button to add some videos to your project. The number of imported videos is not limited. Note that AnimalTA is not doing a copy of these videos, also once a video has been imported in the project if you rename it or change its location you will not be to use it in AnimalTA anymore.

If one or more of the videos you want to import are not in *.avi* format, a window will appear proposing you to convert them. You can then choose which of these videos you want to convert. Only the converted videos will be added to the project. The converted videos may be found in the project folder, in the “converted_vids” subdirectory.

For each imported video, a line will appear in the table presents in the center of the main project (see section: [Main project table](#)). A video can be selected by clicking on its name, in this case,

the rectangle at the left of the video name will appear in green to indicate which video is selected.

2. Import videos from another project

To import an existing AnimalTA project into a new one, you can follow the same steps as you would to add videos. Simply click the "Add video(s)" button and then select the AnimalTA file that you want to import (with a *.ata* extension) instead of video files. By doing this, all the videos that were present in the old file will be imported into the new project along with any associated data, such as parameters and tracking data. This will help you to transfer your existing work quickly and easily to a new project without having to manually add all the videos and associated data again.

3. Suppress video

To suppress a video, click the red cross (✖) under the "Video" label. Note that this action will remove the video from the project and cannot be canceled after validation. All the data associated with this video will also be suppressed (deletion of the tracking parameters and coordinates resulting from tracking). The original file of the video will not be impacted by this action.

4. Duplicate video

Click on the double sheets (📄📄) to duplicate a video. A new video row will appear in your table which the exact same characteristics as the original one, except for its name. This new video is completely independent from the previous one.

5. Concatenate videos

Often, video cameras recording for a long period of time will create several videos of short durations instead of one long video. To add these short videos into one long, you must first import all the parts of the video in the project (see section: [Add videos](#)). Once the importation done, click first on the symbol representing a torn film of the first video (🎞️). The video after which other videos will be concatenated now appears in blue. Click on a second video's name to paste this second video after the first one. The second video will be added after the first one and will disappear from the main project table. You can then repeat this operation as much as you need. If you pressed the "concatenate video" button by accident, you can cancel it by selecting again the same video.

IV. Main project table

In the main window of the project, each imported video is associated with one line of the table presents in the center of the window. A vertical scrollbar at the right of the table allow to scroll through the videos if they are too numerous to fit within the window. The mouse wheel or up and down arrow keys can also be used to scroll. A horizontal scrollbar at the bottom of the table allows to scroll from left to right to see the whole table.

Each row of the table is separated in several sub-elements, described below from left to right. Changing these elements is not mandatory and the program will work even if no changes are done.

1. Video identity

The first section of each row contains three elements:

- A vertical rectangle whose color indicates whether the video is selected (green) or not (red). It is not possible to select more than one video at a time.
- The video name.
- The eye symbol. You can see the first frame of the video by putting your mouse over this symbol.

You can change the name of a video within the project by simply click and write in the video name text field. You cannot have two videos with the same name or with an empty name. AnimalTA will indicate un-accepted names by coloring the video row in red. Changing the name of a video inside the project won't change the file name of the video out of AnimalTA. The new name will be the one used to refer to this video in AnimalTA and in the data created by the program.

2. Frame rate

The fps (frame per second) or frame rate of a video corresponds to the number of images this video is displaying per second. Here reducing the fps means that you will only analyze part of the images of your video. For instance, dividing the fps by two mean that only one out of two images will be analyzed.

By reducing the frame rate, you can greatly improve the speed of the tracking and of the analyses. Decreasing too much the fps will decrease the precision of your measurements and the distances traveled by your subjects will be underestimated. In the opposite, a too high frame rate might result in a confusion between real displacements of the subjects and the noise related to the tracking imprecisions and induce an overestimation of the distances traveled. This last problem can be later corrected by applying a smoothing filter on the coordinates (see section: [Prepare analyses, Smoothing filter](#)).

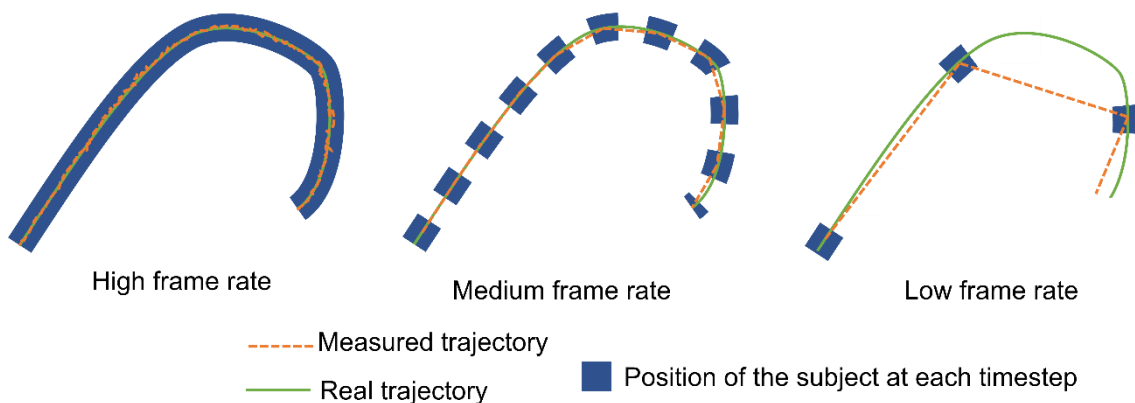


Figure 1: Illustration of the impact of the frame rate on the tracking results.

If you are working with more than one video at a time, you can use the “Apply to other videos” button to use the same frame rate reduction to other videos. Note that the frame rate reduction consists in dividing the number of fps by a given number, also if your videos did not share the same frame rate before the transformation, they will also end up with different frame rates.

If you don't change this option: the original video frame rate will be used.

3. Cropping

Temporal cropping

Many researchers use to record longer movies than needed when recording animal behavior, this method ensure that the whole test will be correctly recorded without missing the first or last seconds of the experiment. This second section allows to remove the useless part of the video *via* cropping.

When you press on the “Crop video” or “Modify cropping” buttons, the main project window will be replaced by a new one. In the center of this new window, you will see the video reader (see section: [video reader](#)) which displays your video.

Under the video reader, you have access to the information about your video cropping both in terms of seconds and frames. To indicate where you want your video to begin, you can either write it in the text field under the “start” label or go to the frame you want as the first frame using the video reader and then press the “Select video beginning” button. The time-bar will be immediately updated, and you will see a blue rectangle indicating where your video will begin (after the cropping). To indicate where you want you video to end you can either click on the time-bar and press the “Select video end” button or indicate the duration you want for your video by writing in the text fields below the “Duration” label.

Spatial cropping

Sometimes, you may also be interested in a reframe of the video. To reframe the video, click on the borders of the displayed image and drag toward the center of the frame. You will see that a red delimitation will appear, only the part inside these delimitations (which will appear brighter than the rest of the image) will be used for the video tracking. You can then repeat the operation with the three other borders. Note that removing this way a significant part of the video will accelerate the tracking process, but removing a very small part of the video might in the opposite slow down the process.

Validate cropping

Once you finished the cropping processes, you can press the “validate” or “validate and go to next video” buttons. The first one will bring you back to the main project table, while the second will open the next video.

You can also move from one video to another using the drop menu in the top left corner of the window. Note that in this case, the changes you made regarding the cropping will not be saved.

If you don't change this option: the whole video will be tracked.

4. Stabilization

This option allows to correct the video in case of camera tremors during the recording. Shortly, this option will look for points of interest in the first frame (after cropping) of the movie. Then, for each frame of the video these same points are identified, and a transformation is applied to ensure that the points of interest remain in the same place through the video. You can ensure that the result of the stabilization is satisfying by clicking the “Check stabilization” button. This option will open the video reader (see section: [video reader](#)) that will show which points of interest have been selected on the frame of reference (first frame). You will also see below the original video with the position of the points of interest and the video after stabilization. If your background is highly homogenous, the program might have difficulties in finding good points of interest. Also, if you see that the stabilized video is not satisfactory, you can remove the problematic points by clicking on them in the reference frame from the video reader. A minimum of two points is mandatory. You can also play with different parameters (in the right part of the windows) to change the way the points of interest are found by the program.

Finally, the “Apply/remove stabilization to multiple videos” button in the main project window allows to generalize this option and its parameters to several videos at the same time.

Note that the process of stabilization is time consuming, also we recommend to not apply a video stabilization on videos that are already stable.

For more technical information, see:

https://docs.opencv.org/3.4/dc/d6b/group_video_track.html#ga473e4b886d0bcc6b65831eb88ed93323

and

<https://learnopencv.com/video-stabilization-using-point-feature-matching-in-opencv/>

If you don't change this option: the original cropped video will be used.

5. Background

To detect the moving targets in the video, AnimalTA uses a process of subtraction: it compares the frame with and without the targets to find the target's position (see Figure 2).

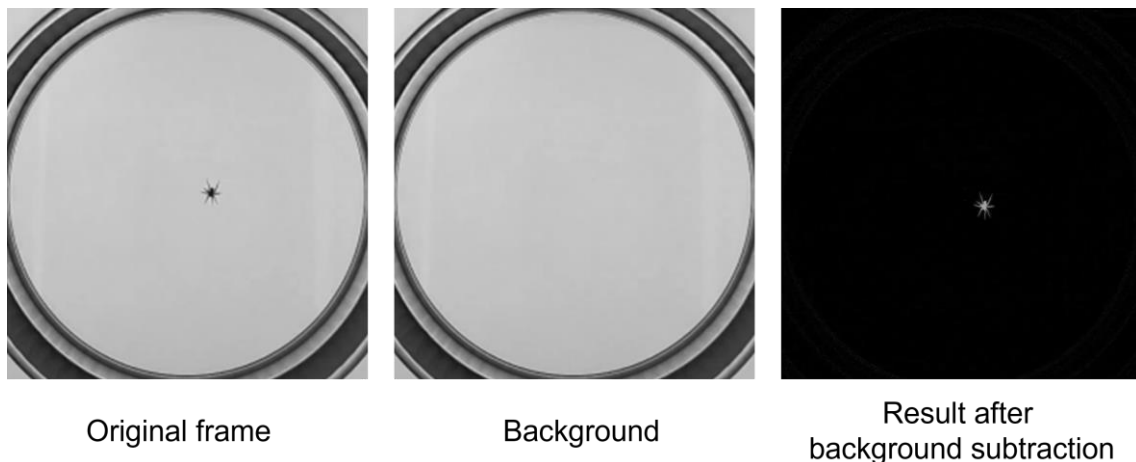


Figure 2: Illustration of the process of background subtraction.

Also, prior to the tracking, it is important to have a background image (grayscale image without the targets). To automatically create this background image, you can press the “Create background button”. To create the background, AnimalTA will take a subsample of the images of the video and calculate the median color value of each pixel of the frame. This methodology will successfully remove all the targets from the image, except if these targets remain at the same place for more than half the duration of the video. In this case, the automatic background process will fail. To ensure that such problem did not arise, you can put your mouse over the eye symbol to see the background proposed by AnimalTA. If some targets are still present, even partially, you will need to correct this background. The correction can be simply done using the “Modify background” option.

This option will open a new window in which you will find the current background of the video, the objective is then to draw over the remaining targets to hide them and to recreate an image empty of targets. To do so, you can select a new color by right clicking in the image, the color under your cursor will be taken (the current color is displayed on the bottom right of the window, above the “validate” button. You can then paint over the targets using left click. The size of the painting area can be adjusted using the mouse wheel. If you are unsatisfied with the automatic background proposed by AnimalTA you can also directly use the first frame of the video by using the button in the bottom right corner.

Keys	Actions
Right click	Change the current color
Left click	Paint
Ctrl + left click on the image	Zoom in
Ctrl + right click on the image	Zoom out
Mouse wheel	Change the size of the tool
Ctrl + Z	Cancel last painting operation

Note that the “Do auto background for multiple videos” button will not reuse the same background for all videos but will recreate a new background of each of the videos.

If you don't change this option: The method used to find the targets will be different and an adaptive thresholding will be used instead of a background subtraction. This second method can be better in case of videos with a lot of changes in luminosity or with important perturbations. With the adaptive threshold method, the program will look for black contrasts on each of the frame to localize the targets.

6. Definition of the arenas

This option allows to draw the delimitations or the arenas in which the targets may be found. Drawing the arenas ensures you that the targets will not be confused between arenas during the tracking: a target inside an arena cannot leave this arena or appear in another one. It also limits the risks that some elements outside the arena will be confused with the targets. Finally, the defined areas will be later used in the analyses to calculate for instance the proportion of the arena explored by the target, or the time spend close to the border of the arena.

Once you pressed the “Define arenas” button, a new window will appear. In this window you will see the image of your background (first image of the video if the background has not been defined) and you will be able to draw the different arenas. To draw an arena, simply click on the

background where you want to place a corner of the first arena. At any moment, you can indicate the kind of shape your arena should be (Ellipse/Rectangle or Polygon). Once the corners are placed, you can move them by click and drag, or suppress them using the “supr” key. Once your first arena is correctly defined, you can press the “enter” key to validate it and begin the drawing of a new one. You can modify a precedent shape simply by clicking on one of its points. Arenas can also be moved, rotated, or resized as described in the table below. One area must be completely separated from the others, also if two shapes are overlapping, the combination of these two shapes will be considered as a unique area. You can use the “View” button to see the result. If the arenas have the same location in several videos, you can use the “Apply these arenas to other videos button”. In case of doubt the number of arenas is displayed below the information panel.

Keys	Actions
Left click	Add a corner
Left click on a point	Select/move the corner
Left click and drag on a border	Move the arena
Suppr	Delete the selected corner
Enter	Validate a shape
<Shift> + Left click and drag on a border	Duplicate the arena
<Shift> + Right click and drag on a border	Resize the arena
<R> + Right click and drag on a border	Rotate the arena
Ctrl + left click on the image	Zoom in
Ctrl + right click on the image	Zoom out


If you don't change this option: The program will consider that there is only one arena occupying the whole image.

7. Definition of the scale

Defining the scale of your recording will allow you to obtain data expressed in international units instead of pixels. Another advantage of defining the scale of your video is that it will allow to uniformize the tracking parameters between videos which are not recorded at the same scale (e.g. Figure 3).


After clicking the “Define scale” button, a new window will open in which you can see the first frame of the video. Click on the frame to place two points separated by a known distance. Then, indicate in the bottom right files the real distance between these two points and which unit you are using. If you already know the scale of your video, you can directly type in the “Ratio” text field.

You can use the “Apply this scale to other video” if several videos have been recorded at the same scale (same distance between the camera and the arenas, same quality of recordings, same zoom).

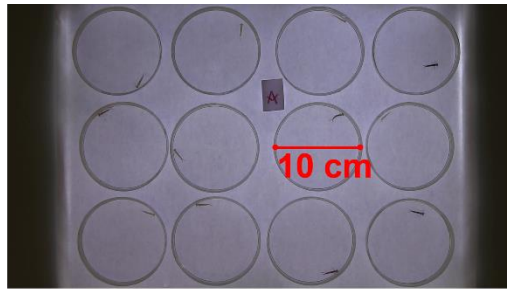

Video 1:

Mean target length: 35 pixels / 1.5 cm

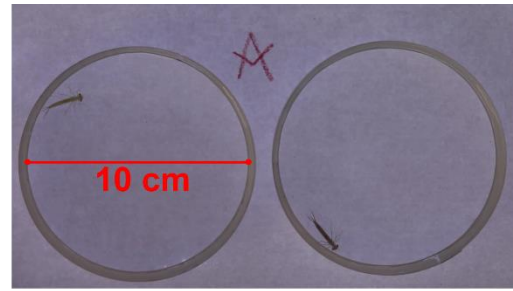

Video 2:

Mean target length: 135 pixels / 1.5 cm

Figure 3: Illustration of the importance of the scaling: when the recording is not done from the same point of view, the values in pixels cannot be compared between videos but the distance in international units can be.

If you don't change this option: All the later measurements will be done in pixels.

V. Video reader

The video reader will be found in different situations using AnimalTA. It will show you the video you are working with (Fig 4, A). Under the video, you can see a time-bar (Fig4, B) that will show you the duration of the video and which indicate which frame of the video is visible on the screen (indicated by a black rectangle, Fig4 C). The red part of the time-bar is the part of the video you are working with (i.e. video after cropping, see section: [Main project table, Cropping](#); Fig 4, D), it is delimited by a blue and a red rectangles who indicate respectively the beginning and the end of the working part of the video. You can click directly on this time-bar to go to a particular moment of the video.

Four buttons are available under the time-bar:

Play/Stop buttons: play/stop the video

Go to begin/end buttons: Make you jump at the beginning/end of the video

A slider at the right of these buttons allows you to increase/decrease the speed at which the video is displayed while playing. A speed of 0 corresponds to the original video speed.

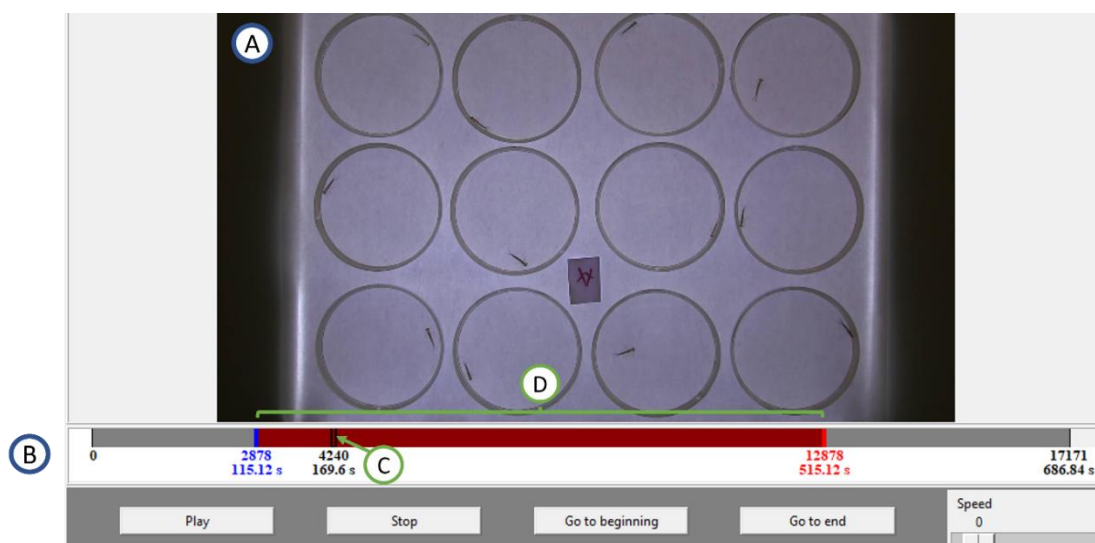


Figure 4: Illustration of the video reader.

Shortcuts:

Keys	Actions
Ctrl + left click on the image	Zoom in
Ctrl + right click on the image	Zoom out
Spacebar	Play/Stop the video
Left arrow	Move backward of one frame
Right arrow	Move forward of one frame

VI. Tracking preparation

Once you finished with the preparation of the videos (see section [Main project table](#)), you will be able to define the parameters of the tracking. To that aim, select the video you are interested in and press the “Prepare tracking” button in the bottom left of the main project window. This will direct you toward a new window in which you will see the video reader (see section [video reader](#)) and several parameters in the right of the window. Each parameter corresponds to a transformation made to the video that will allow to detect the targets in the image. The transformations will be applied from top to bottom, also if you select the 4th checkbox (Erosion) it means that the image displayed in the video reader is the image that went through the 3 other transformations (Original video importation, Background subtraction and Thresholding).

From the main project window, you can see the tracking preparation’s parameters of a selected video in the bottom right of the window. You can also apply these parameters to other videos using the “Apply these parameters to other videos” button.

Original video

Show the video cropped and stabilized if this option was selected.

Gray image

Show the video in grayscale mode. If activated, the “correct brightness” button induces a correction of the brightness. To that aim, the brightness values of the current frame are modified according to the distribution of a reference image brightness values. The reference image is the background, if there is no background, the reference image is the first frame of the video.

Background

If a background has been defined (see section [Main project table, Background](#)), this option shows the video after the background subtraction. Note that the process is not a simple background subtraction process but a combination of two subtractions ($\text{sub}(\text{background}, \text{frame}) + \text{sub}(\text{frame}, \text{background})$) to detect both light and dark targets.

Threshold

If a background has been defined, this parameter transforms the image resulting from the background process to create a binary image (only composed of white and black pixels, no grey anymore). The threshold parameter determines the lightness above which a grey pixel will be converted to a white one and below which it will be converted to a black one (Figure 5). For a good quality tracking, the result (visible in the video reader) must be the closest possible of the targets in white and the rest of the image in black.

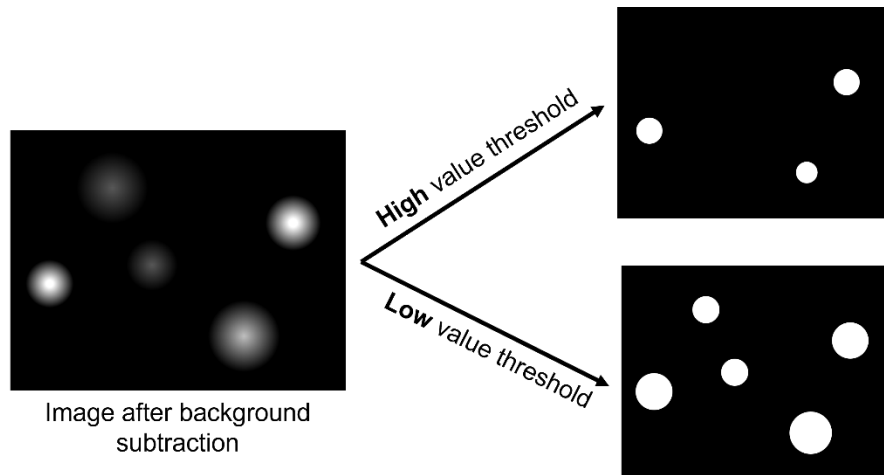


Figure 5: Illustration of the process of thresholding.

If there was no background defined for the selected video, this process will be replaced by an adaptive thresholding: the program will identify spots that are darker than their close environment. The parameter can be changed to change the range of considered environment (Figure 6).

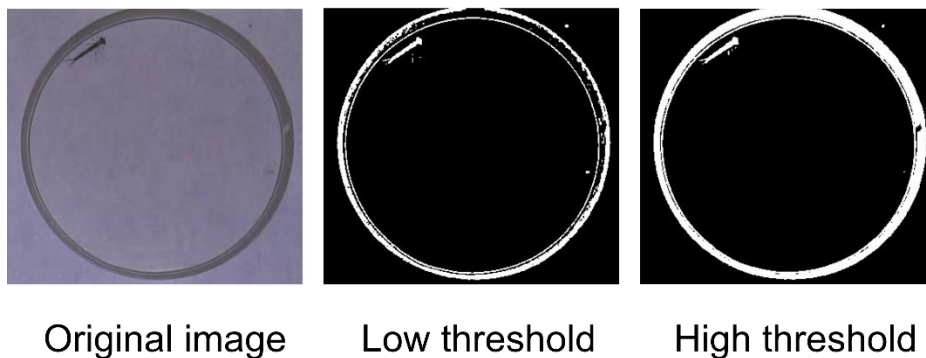


Figure 6: Illustration of the process of adaptive thresholding.

Erosion/Dilation

This allows to increase/decrease the size of the white areas resulting from the threshold process. While erode will remove white pixels from the borders of the white shapes, dilate will add some. These parameters allow to remove potential noise from the background and improve the visibility of the targets (see Figure 7).

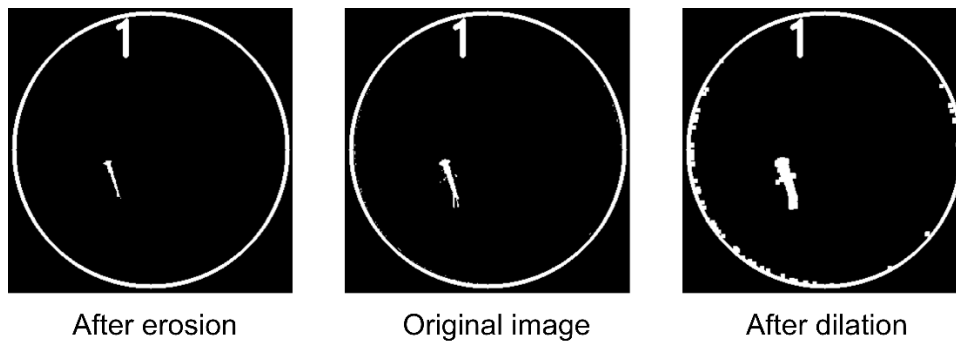


Figure 7: Illustration of the processes of erosion and dilation.

Area selection

This parameter allows to determine the minimum and maximum area of your targets. All white shapes that are smaller/bigger than the indicated value will be removed from the image (Figure 8). This way you will limit the risk that another element is considered as a target.

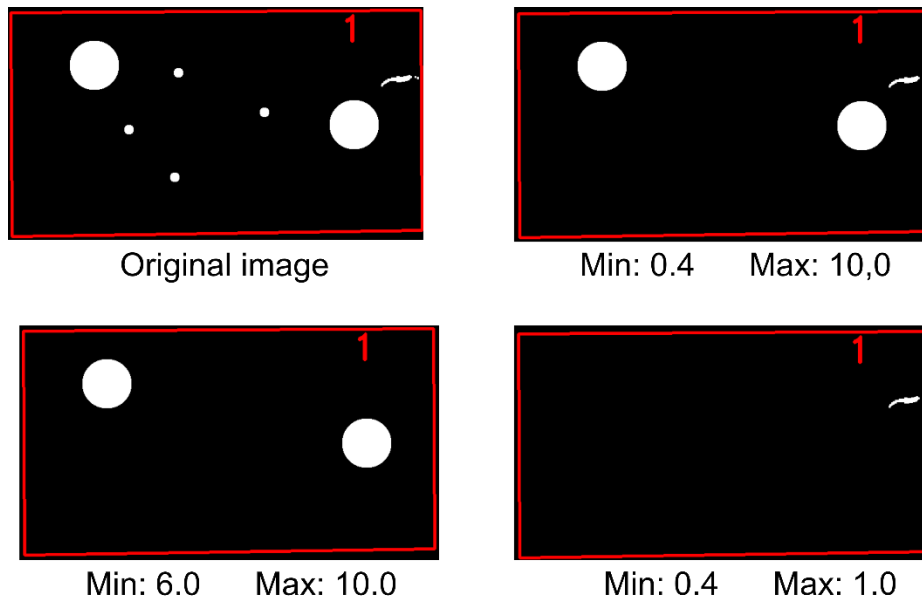


Figure 8: Illustration of the process of area selection. Under each image is the value set for minimum and maximum areas.

Distance threshold

This parameter may be used to indicate the maximum distance a target can move between two frames. During the tracking process, the distance between the target at a time t and the potentials targets at $t+1$ is calculated, and all the potential targets whose distance from the target position at time t is higher than the threshold distance will be excluded. The best value for this parameter is not easy to find as it will depend on a lot of parameters like the maximum moving speed of the targets or the frame rate of the video.

To facilitate this process, after clicking the “Distance threshold” checkbox, you will see in the video reader light blue points: they indicate where the potential targets are located. You will also see dark blue disks whose center is the position of the potential targets at the previous image and whose radius is your distance threshold. Assuming that a target was correctly located at time $t-1$, when at t the target light point is inside the corresponding dark circle, it means that AnimalTA will consider this light point as a potential position of the target for time t (Figure 9). We advise to look for moments of the video where the targets are moving fast and to ensure the distance threshold is not too low. It will always be better to use a too high than a too low distance threshold.

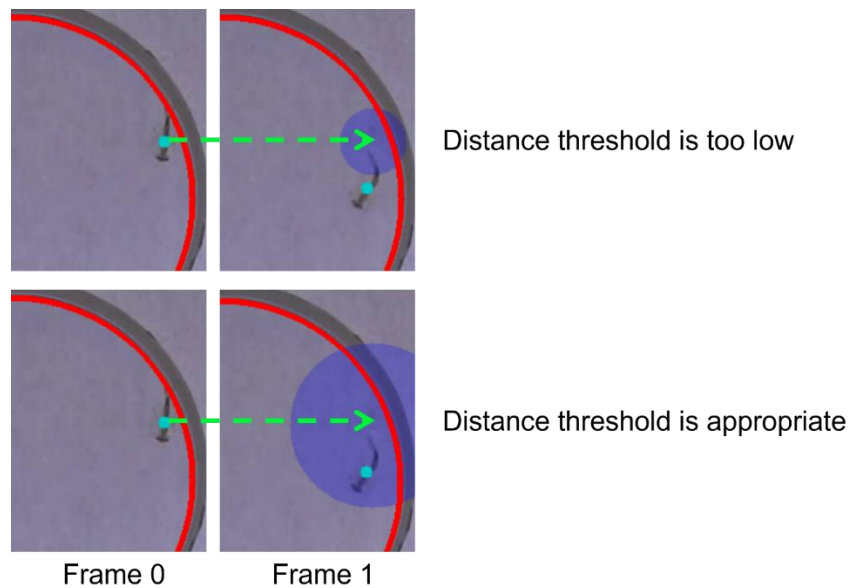


Figure 9: Illustration of the distance threshold.

Number of individuals per arena

According to your experimental design, you might have more than one target per arena, in that case you must indicate how much targets are to be found in each area. If this number is the same for all the arenas, you can just change the slider to adapt it to your needs. On the image, each arena is drawn in red, and the number of targets expected to be found appear as a red number close to the arena border.

If the number of targets per arena is variable, press the “P” button at the right of the “Number of individuals per arena” section to manually change the number of individuals for each arena.

Note that AnimalTA does not allow to track the identity of the different targets within a same area. For such kind of analyses, please see: <https://www.idtracker.es/> from Pérez-Escudero et al., 2014.

A. Pérez-Escudero, J. Vicente-Page, R.C. Hinz, S. Arganda, G.G. de Polavieja, *Nature Methods* 11(7):743-748 (2014)

VII. Tracking process

Once that at least one of the videos is ready for the tracking, you will be able to click the “Begin tracking” button and to select the videos you want to track. According to the number of videos, their length, their quality, and the options your selected, this process duration can be highly variable.

VIII. Visualization and correction of the trackings

1. Visualization

You can check the accuracy of the trackings using the “View and correct track” button in the bottom right part of the main window. A new window will then appear with the video reader on the left on which you can see the results of the tracking. The coordinates of the targets are marked by a colored point and the tails behind these points represents the previous positions of the targets. At the same time, you can see the coordinates in a table at the right of the window (Figure 10). The length of the tail may be increased using a slider which represent the tail length in seconds (Figure 10, A). You can also make appear all the trajectories recorded during the video by pressing the “see all trajectories” button (Figure 10, E).

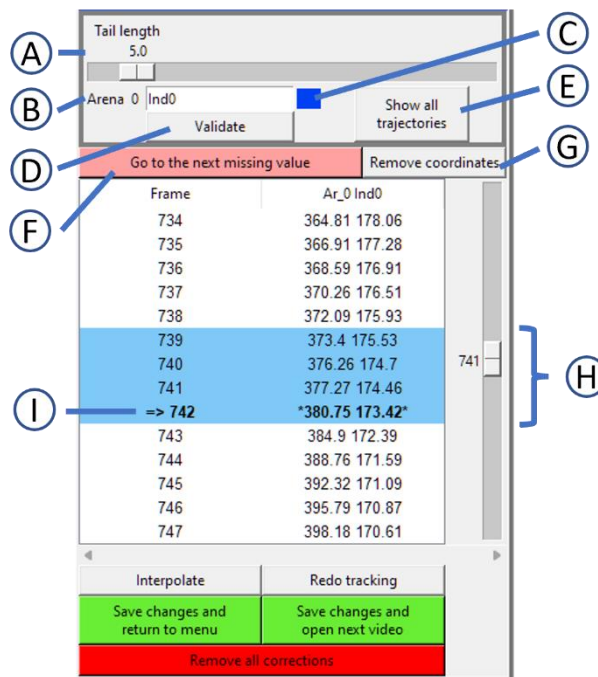


Figure 10: Correction tracking panel. A) length of the tail displayed, B) Name of the target, C) Color of the target, D) Validate the name of the target, E) Draw all trajectories on the frame, F) Move the video up to the next missing value, G) Replace coordinates by NA values, H) Selected frames, I) Current frame.

In the table at the right of the window, you can see the coordinates corresponding to what is displayed in the video reader. The frame displayed is indicated by “=>” preceding its number (Figure 10, I). Each column of the table corresponds to one target, each cell contains the X and Y coordinates of these targets. In addition to using the features of the video reader (see section [Video Reader](#)), you can set the video reader to a particular frame by clicking on a row from the table.

2. Set the identities of the targets

By default, the targets are identified with an unrepeated within area label, going from “Ind0” to “IndN”, N+1 being the number of targets expected inside the arena (see section [Tracking preparation, Number of individuals per arena](#)). When you select an individual (by pressing on it in the image or by pressing on its column in the trajectories table), its identity will be displayed in the text field on the top of the trajectories table (Figure 10, B). You can change the label of the selected individual by changing the text in the corresponding text field and then press the validate button (Figure 10, D). The new label will later be used when performing data analyses. The color used to display the position and trajectory of the selected target can also be modified by clicking the colored rectangle at the right of the text field (Figure 10, C). A color chart will then appear, just click on the color you want to use, and this new color will be associated with the target.

3. Corrections

There are three different ways of correcting tracking errors that must be used according to the duration of the according error.

Find missing values

If the targets are lost during a part of the tracking, their coordinates will appear as “NA” values. To correct rapidly the video, you can press the “Go to the next missing value” button (Figure 10, F), which will transport you directly to the next problematic frame. If your dataset has no missing value, this button will appear green and cannot be clicked.

Punctual error

To correct a punctual error which happened during only one frame, you can directly change the position of the points on the video reader by click and drag it, the table will be updated accordingly.

If there are no points at all (i.e. NA value), you can directly click right on the frame where you want to add the point, and the coordinates of the target will be added to the table. To facilitate these manual corrections, after each right click the video will be moved forward by one frame, allowing to use only right clicks to make frame by frame corrections.

Interpolation

A single target. If the target has been lost during a short duration and that its trajectory during this period can be resumed by a straight line, you can interpolate the coordinates between two frames (see Figure 11). To do so, you must first select the concerned target. You can either click on it on the video reader or click on one corresponding cell of the table. The selected individual may be identified on the video as the point will be highlighted, on the table, the corresponding coordinates will appear between “ * ”. You must then select the frames on which you want to apply the interpolation (Figure 10, H). To that aim, you can either press <Shift> and click on another frame of the table or play the video maintaining the <Shift> key pressed. Once you made this selection, press the “Interpolate” button to apply the changes to the trajectories.

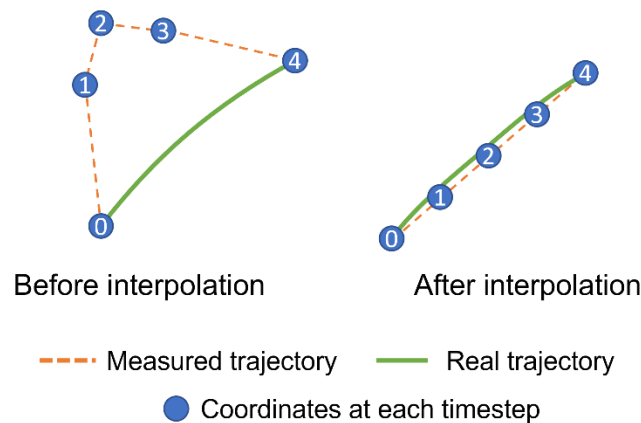


Figure 11: Illustration of the process of interpolation. In this case, an interpolation is realized between frame 0 and 4.

All targets. If you want to fill missing values of all targets for the whole video, you can press the “Interpolate” button when less than two frames are selected in the data frame (see Fig 10, H; section *Interpolation*). This will open a warning window, just click “yes” and the data will be immediately updated.

Redo part of the tracking

If you need to correct some bigger issues, you can simply redo a part of the tracking with different parameters. To do so, select the frames you want to change (see section *Interpolation*) and press the “Redo tracking button”. This will open a new window in which you will be able to choose the pre-tracking parameters for this portion of video (see sections [Main project table](#) and [Tracking preparations](#)). After changing the parameters, press the “Redo the track” button to see the result. Once you are satisfied with the result, press the “Validate these changes” button to apply them within the rest of the video.

Invert coordinates

If the identity of two targets has been inverted, you can swap it back by selecting one of the two targets and then click on the second one (either on the image or in the table) while holding the “Shift” key. The coordinates of these two targets will be inverted from the current frame until the end of the video.

Delete existing values

If you want to replace some existing coordinates by NA values, you can select them (Figure 10, H; section *Interpolation*) and press the “Remove coordinates” button (Figure 10, G).

Validate corrections

Once you pressed one of the two validate button, a new coordinates file will be created inside the project folder (see section [Generalities and project management](#)) under the “corrected coordinates” subdirectory. You can also remove all the corrections you did by pressing the “Remove all corrections” button.

IX. Prepare analyses

From the main Project window, you can press the “Prepare analyses” button in the bottom right of the window after selecting a video.

1. Smoothing filter

You can smooth the trajectories of your targets by checking the “Smooth trajectories” checkbox. This parameter applies a Savitzky-Golay filter (function `savgol_filter` from the `scipy` python package) on the coordinates. The application of this filter allows to facilitate the identification of the movement phase and to correct the noises associated with the imprecision of the tracking (Figure 12, see also section [Main project table, Frame rate](#)). The parameters of the smoothing filter can be controlled by pressing the “P” button at the right of the checkbox. To compare the between video’s information, you might want to apply the same smoothing filter to all of them at the same time. To do so, you can use the “Apply smoothing to other video(s)” button to apply the smoothing parameters defined for this video to other ones. If no smoothing was applied at this moment, the smoothing filter will be removed in the other designated videos.

For more information, see:

https://docs.scipy.org/doc/scipy/reference/generated/scipy.signal.savgol_filter.html

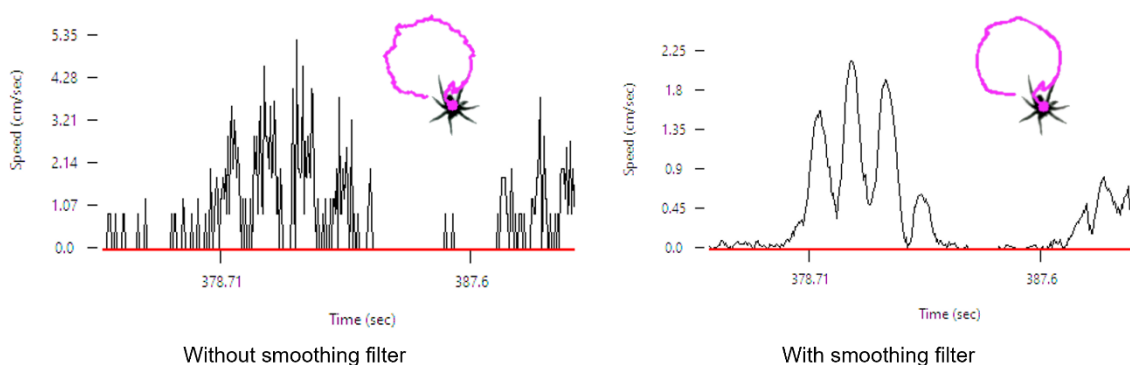


Figure 12: Illustration of two trajectories before and after the application of the smoothing filter. The graphs show the measured speed of the individual as a function of time. We can see that the periods of movement are easier to identify after the application of the filter.

2. Perspective correction

This option allows to correct for camera deformation and/or the effect of perspective. If you press this button, a new window will appear. In this window, you will see a frame of reference (Figure 13, A), a schematic representation of your environment (Figure 13, B), and the result of the deformation (Figure 13, C). You must first click on the reference frame (A) and define at least 4 points of reference. These points will appear in the schematic representation. Move the

points by click and drop inside the schematic representation so that all distances and angles between points are correct. You can delete an existing point by pressing the <Delete> key of your keyboard after selecting it. After validating, you can apply the same movement threshold to all videos using the “Apply the exploration parameters to other video(s)” button. Note that the transformations of the images and trajectories are done after the tracking process, so they won’t impact the tracking quality.

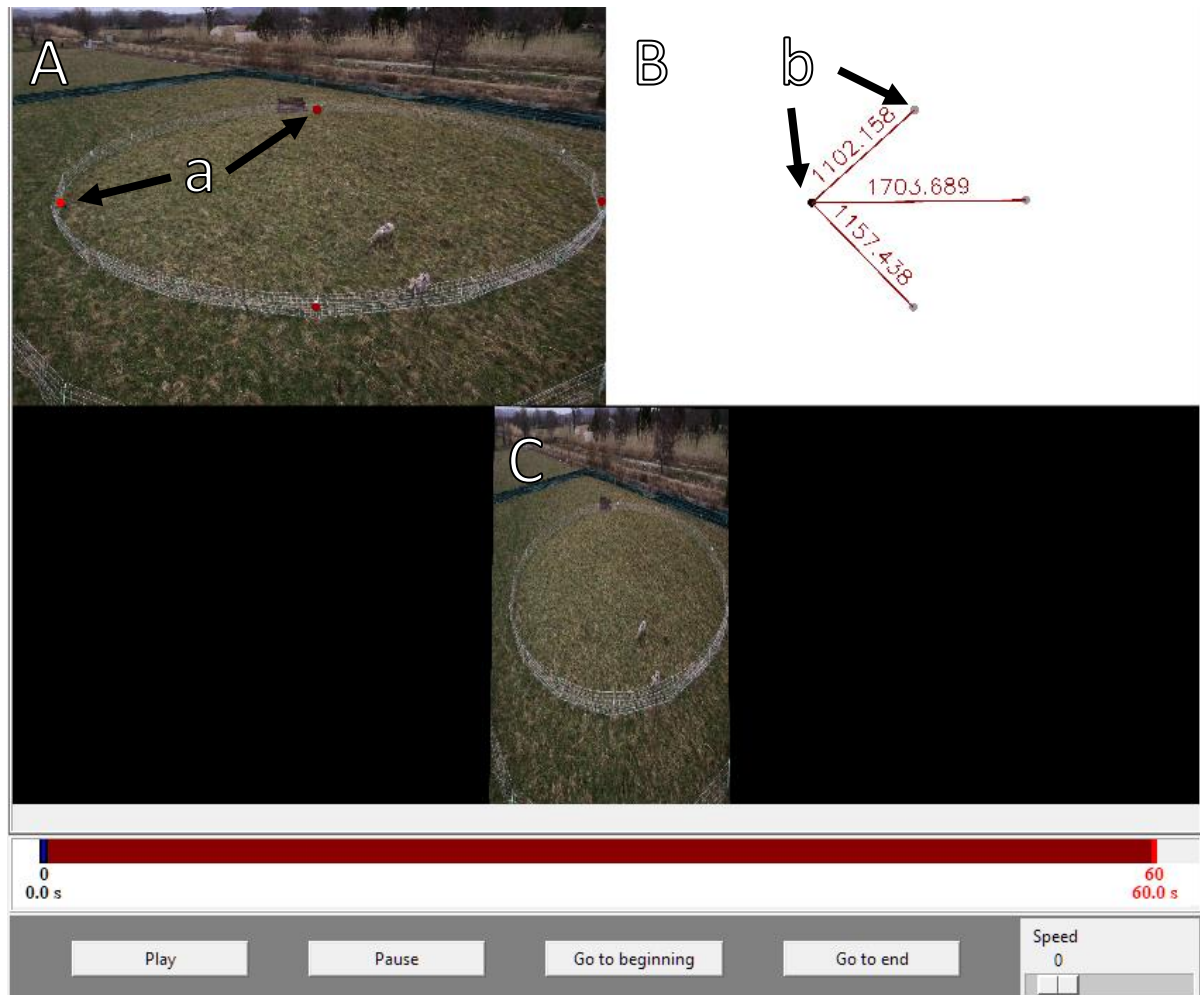


Figure 13: Panel used to correct the effect of perspective. We can see in the reference frame (A) the points of reference (a). A schematic representation of the scene (B) in which the points are also drawn (b). Finally, the result can be seen in the last image (C).

3. Analyses of individual trajectories

If you click the corresponding checkbox, you will see on the video reader the movement speed of each target (except for the first frame). The speed value will appear in green/red if the target is considered as moving/stopped respectively. By default, a target will be considered as moving if the distance traveled between two frames is higher than 0. However, this method does not consider the tracking imprecisions which normally results in little trembling of the target.

To change the threshold of movement above which you want to consider that a target is moving, you can press the “P” button. This will open a new window in which you can see a graph of the speed of the selected target as a function of the time. You can press the drop menu on the top

of the window to change the selected target, the selected target will also be highlighted on the video reader. A blue vertical line on the graph indicates the currently displayed frame in the video reader, you can move this line by simply clicking somewhere on the graph; the video reader will be immediately updated. A red horizontal line indicates the movement threshold: the speed above which an individual will be considered as moving. You can change this value by click and drag the red line.

Finally, in the bottom right corner of the window, you will have a summary of the measured data for the selected target:

- Proportion of time lost: Count the proportion of frames for which the target was not detected (appears as "NA NA" in the data tables).
- Average speed: The displacement speed between each frame will be calculated and averaged for the whole video.
- Average speed when moving: The displacement speed between each frame where the target is considered as moving (i.e. speed above the red movement threshold) will be calculated and averaged for the whole video.
- Proportion of time moving: Proportion of time the target is considered as moving.
- Total distance traveled: Sum the total distance traveled by the target.
- Total distance traveled when moving: Sum the total distance traveled by the target when we only consider the frame during which the target is considered as moving.

After validating, you can apply the same movement threshold to all videos using the "Apply the exploration parameters to other video(s)" button.

4. Spatial repartition

This section of AnimalTA allows to define some points or areas of interest. Click the corresponding checkbox to see the elements that are already defined (non by default). Click the "P" button to add some new elements. This will open a new window in which you can see the current frame of the video with the selected area highlighted. Use the drop menu in the top of the window to change the selected area. You can then press the "Add element" button (top left of the window) to add new elements.

Elements manipulation

The elements placed on the image can be easily moved by click and drag the associated points (except for borders, see below). When created, each element will appear in the right part of the window where you can see its name (by default the kind of element and an identity number) and the measurements associated with this element. You can rename an element by simply writing in its name text field. Note that you cannot have two elements with the same name. You can also suppress an element by pressing the "Suppr" button at the right of its name.

Types of elements

- Point

A point can be either freely placed (just click on the frame where you want to place it) or placed in the center of the area.

A point is associated with three different measurements: i) the average distance to this point, ii) the latency to approach at less than a given distance to this point (to be defined by the user by

simply writing in the correspondent text field), iii) the proportion of time the target spent inside this at less than this distance from the point.

- Segment

A segment is defined by two points. These two points can be placed freely (click anywhere on the frame) or be associated with a border. If you choose the second option, you must click on the border/corner of the area where you want to place the first point of your segment. If you select a corner, the point will be immediately placed on the frame. If you choose a border, a small window will appear that will ask you what position you want this point to take on the border. Then repeat the same operation for the second point.

Three measurements are associated with a segment. i) The mean distance to the segment which is the average value measured for all frames as the shortest distance between the target and any point of the segment. ii) The number of times the target crossed the segment. iii) The latency before crossing the segment for the first time.

- Borders (All)

This will immediately select all the borders of the area (i.e. the whole contour of the area) and add it as a new element.

Two measurements are associated with this. i) The mean distance to the border is the averaged shortest distance between the target and the border. ii) The proportion of time the target spent at less than a given distance from the border (to be defined by the user by simply writing in the correspondent text field).

- Borders (selection)

You can also select specific borders. In this case, you must indicate on the frame which border(s) you want to consider. To do so, click on the border(s) you want to consider (their color will change from purple to red if they are selected). Once you selected all the borders you want, press the "Return" key to validate and the borders will appear in the list of elements.

Three measurements are associated with this. i) The mean distance to the border is the averaged shortest distance between the target and the selected border(s). ii) The proportion of time the target spent at less than a given distance from the border (to be defined by the user by simply writing in the correspondent text field). iii) The latency to approach the border(s) at less than the indicated distance.

- Area

To add an elliptical, circular or polygonal area, you must click on the frame where you want the corners of your area to be placed. Once all the corners are placed, press the "Return" key to validate the shape, which will appear in the list of elements.

Two measurements are associated with shapes: i) the proportion of time the target spent inside this area and ii) the latency to enter this area.

Expend elements to other arenas

Once you defined at least one element, you can duplicate this element and apply it in other arenas. To that aim, click the "Expend to other arena(s)" button. This will open a new window with two lists on it. The left list is the list of elements present in your arena; from this list you can select the elements you want to duplicate in other arenas. The second list is the list of arenas (all videos confounded). If you want to visualize an arena, just put your cursor over its name in

the list, an image of the video with the arena drawn in blue will appear. Select in this list the arenas toward which you want to duplicate the elements. Simply click the “Validate” button to apply these duplications. It is important to note that except for the “All borders” elements the other elements will be duplicated only if the arena of origin and of reception are of the same shape (same number of corners, square/circle, ...). If the shapes are too different or do not share the same number of corners, the duplication of elements will not be applied. Moreover, two elements in the same arena cannot have the same name. Also, if you try to add a new element in an arena which already possessed an element with the same name, the old element will be replaced by the new one.

5. Interindividual distance

The interindividual distance is simply measured as the sum of all distances between all the targets inside an arena. By clicking this checkbox, you will see on the video reader the calculated distances for the current frame.

You can also define a threshold distance by using the “P” button at the right of the checkbox. This will open a new window in which you can set a threshold either by moving a slider or by directly writing in a text field. If the distance between two targets at a given timestep is lower than this threshold, these targets will be considered as “neighbors” when you will be running the analyses (see section [Run analyses, Results of interindividual distance](#)).

6. Exploration

This option allows to calculate the proportion of the arena explored by each target. To change the parameters associated with the measurement of exploration, you can choose between them clicking the “P” button at the right of the checkbox. This will open a new window in which you can select a particular individual (you can change the current individual using the drop menu at the top of the window). AnimalTA proposes two different methods to measure the exploration. After defining the method of exploration, you can apply it to other videos using the “Apply the exploration parameters to other video(s)” button.

Individual centered

This method allows a fine measurement of the explored area. With this method, we consider that a target explores a circular area around it at each frame. Each pixel that falls within this circular area at least once during the video is considered as visited while the others are considered as unvisited. The measure of exploration returned is simply the proportion of visited pixels within the arena. You can change the surface of the circular area around the target using the slider or by directly writing in the corresponding text field.

Meshes

This method is more traditional as it has been primarily invented to analyze videos by hand. It consists in drawing a mesh on the top of the arena and counting the number of visited cells of the mesh. We propose here two types of mesh: a rectangular or a circular one. In both cases, the surface of each cell of the mesh can be determined by the user. Regarding the circular mesh, a supplementary parameter can be defined (aspect parameter) which determine the shape of the mesh cells. To method used to draw the circular mesh are based on the work of Beckers and Beckers: Beckers, B., & Beckers, P. (2012). A general rule for disk and hemisphere partition into equal-area cells. *Computational Geometry*, 45(7), 275-283.

X. Run analyses

After you prepared the analyses, you can run them using the “Run analyses” button in the main project window. The results will be summarized as .csv files in the project folder, in the “Results” subfolder. In this subfolder, you will find two .csv files.

Results per individual

The file “Results_by_ind” is constituted of 16 columns and of one row per tracked target.

Column name	Content	Supplementary information
Video	The file name of the video	
Arena	The ID number of the arena in which the target was (the arenas are numerated from left to right and from top to bottom).	Section Main project table, Definition of the arenas
Individual	The ID of the target, if not modified by the user: from Ind0 to IndN, unique for each arena. The user can also personalize the target’s ID (see section Visualization and correction of trackings, set the identities of the targets).	If two targets have been confounded during the tracking process, you can restore the target’s identity following the steps described in section Visualization and correction of trackings, Corrections, Invert coordinates
Prop_time_lost	Calculate the proportion of frames for which the target’s coordinates were unknown.	In case of tracking errors or lost targets, you can correct the trajectories (section Visualization and correction of trackings, Corrections).
Smoothing_filter_window	The window parameter set for the smoothing filter. If the value is 0, it means that no smoothing filter was applied.	See section Prepare analyses, Smoothing filter .
Smoothing_Polyorder	The polyorder parameter set for the smoothing filter. If the value is 0, it means that no smoothing filter was applied.	

Moving_threshold	The speed threshold set by the user above which an individual is considered in movement.	See section Prepare analyses, Analyses of individual trajectories.
Prop_time_moving	Proportion of frames the target is moving faster than the moving threshold.	
Average_Speed	The averaged displacement speed of the target between each frame.	
Average_Speed_Moving	Similar to “Average_Speed” value but considering only the frames where the target was moving faster than the moving threshold.	
Traveled_Dist	Total distance traveled by the target.	
Meander	Estimation of the meander of the trajectories. This value is the average of the change in direction (turning angle) divided by the distance traveled for each frame.	
Traveled_Dist_Moving	Similar to “Traveled_Dist” value but considering only the frames where the target was moving faster than the moving threshold.	
Meander_moving	Similar to “Meander” value but considering only the frames where the target was moving faster than the moving threshold.	
Exploration_value	The exploration score of the target. The nature of this value will change according to the method of exploration measurement choose by the user.	See section Prepare analyses, Exploration.
Exploration_method	The method of exploration measurement set by the user.	
Exploration_area	The area used as an exploration unit. If the method used is squared or circular mesh, this value is the surface of each cell of the mesh. If the method is the individual centered one, this value indicates the surface a target cover around it at each moment.	
Exploration_aspect_param	The aspect parameter used in case of circular mesh method. This value will be NA if another method was used.	See section Prepare analyses, Exploration, Meshes.
Mean_nb_neighbours	Calculate the number of neighbors the target has for every frame. These values are then averaged for the whole video. The distance threshold to consider two targets as neighbors is defined by the user.	See section Prepare analyses, Interindividual distance.

Prop_time_with_at_least one_neighbour	Calculate the proportion of frames for which the target has at least one neighbor.	
Mean_shortest_dist_neighbour	For each frame, the distance between the target and the closest neighbor is calculated. These values are then averaged for the whole video.	
Mean_sum_distances to_neighbours	For each frame, the distances between the target and all the other targets (within the same arena) are summed. These sums are then averaged for the whole video.	

Results of interindividual distances

This file is a summary of the inter-individual distances, it contains one row per arena which contained more than one target. The data found here are the mean inter-individual distance (see section [Prepare analyses, interindividual distance](#)). It also computed the smallest interindividual distance measured and the biggest one. Note that the identity of different targets within a same area may be confused when these targets cross each other's, in such situation, the user can manually correct such mistakes (see section [Visualization and correction of the trackings, Corrections, Invert coordinates](#)).

In the "InterInd" folder, you will find one folder per video which contains detailed results for each of the arenas (with at least two targets). For these results, a contact event is defined as continuous period during which two targets are neighbors (separated by a distance lower than the threshold defined by the user; see section [Prepare analyses, Interindividual distance](#)). The results are separated in four .csv files:

- Contact occurrences: a table of contingency recording the number of contact events between the different targets.
- Contact events: a table listing all contacts (see the above definition of contact). Each row of this table corresponds to one contact event and summarize which were the two targets involved (columns "Ind_P1" and "Ind_P2"), the duration of the contact event (column "Duration"), and the moment at which the contact began (column "Beginning").
- Distances: a cross table in which the distances between targets are summarized (averaged for the whole video).
- PropTime: a cross table with the proportions of frames the two targets spent as neighbors.

Spatial

In the spatial folder, you will find one file per element of interest (see section [Prepare analyses, Spatial repartition](#)). The structure of these table is the same as the "Results per individual" ones, the measured variables are described in section [Prepare analyses, Spatial repartition](#).

Another file called “General” can also be found. This file summarize for each element (except elements of the type lines) the minimum, maximum and average number of targets present inside the element during the video.

Detailed data

If you are interested in more detailed analysis, you will find in this directory the details of each target situation per frame. Files will be organized in subfolders whose names are the ones of the videos. The files names refer to the name of the arena in which target is situated and then the name of the target itself. These .csv files contain one row per frame and at least 7 columns. The columns are the following:

Time: from 0 to the total duration of the video, time is expressed in seconds.

X/X_smoothed and Y/Y_Smoothed: The coordinates of the target after transformation (after being scaled and after [smooth-filtering](#) if applicable). The unit of these values is the one defined by the user for the [scale of the video](#). If no unit was defined, the values are expressed in pixels.

Speed: How fast the target moved between the last frame and the current one (in [unit defined by user](#) per seconds).

Moving: Whether the target is considered as moving or not. A target is considered as moving if its speed is higher than the [movement threshold](#) defined by the user (default value is 0).

Distance: The distance traveled by the target between the last and the current frame. The unit is the one defined by the user for the [scale of the video](#). If no unit was defined, this value is expressed in pixels.

Dist_to_X: In the case of arenas with more than one individual, this column gives the distance (in [unit defined by user](#)) between the considered target and the target X.

Dist_to_Y: In the case of arenas with [elements of interest](#) defined by user, this column gives the shortest distance (in [unit defined by user](#)) between the considered target and the element Y. If the element is a shape, a negative value means that the target is inside the element of interest.

XI. Export a video

As AnimalTA indirectly allows video-editing, you can export the videos you used during the program using the “Export the video” button in the main project table. This will open a new window in which you can select at which step of the process you want your video to be saved. You can choose the length of the tail to be displayed, [smoothed](#) or raw trajectories and with [perspective correction](#) or not. Finally, the [target’s identities](#) can be displayed.